

# Northwest

Life Science Specialties, LLC

Premier Products for Superior Life Science Research

## *NWLS<sup>TM</sup>*

# *Glutathione Assay*

---

**Product NWK-GSH01-1 (96 Well Version)**  
***For Research Use Only***



Simple assay kit for quantitative measurement of glutathione in biological samples such as whole blood, tissue homogenates and cell lysates.


## Table of Contents

---

<b>Section</b>	<b>Page</b>
Introduction	3
Intended Use	3
Test Principle	3
General Specifications	4
Kit Contents	4
Required Materials Not Provided	4
Required Instrumentation	4
Warnings, Limitations, Precautions	5
Storage Instructions	5
Assay Preparation	5
Reagent Preparation	6
Sample Handling/Preparation	6
Assay Protocol	8
Data Analysis	8
Performance Details	11
Things To Note	12
References	13
Statement of Limited Warranty	13
Procedure Checklist	14
End-User Notes	15

**Introduction:**

Glutathione (GSH,  $\gamma$ -glutamylcysteinylglycine), the primary non-protein sulfhydryl in aerobic organisms is synthesized in most cells. The ubiquitous tripeptide is formed by the ATP dependent condensation of glutamic acid and cysteine, catalyzed by  $\gamma$ -glutamylcysteinyl synthetase. Glycine is then added by glutathione synthetase to form GSH.

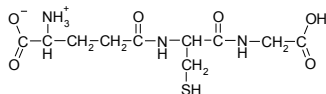

**Figure 1.**  
**Structure of GSH**


In addition to donating an electron during the reduction of hydroperoxides to the respective alcohols (or water in the case of hydrogen peroxide), GSH functions as a co-substrate in the metabolism of xenobiotics catalyzed by glutathione S-transferases. It is also a co-factor for several metabolic enzymes, is involved in intracellular transport, functions as an antioxidant and radioprotectant and facilitates protein folding and degradation.

**Intended Use:**

The NWLSS™ Glutathione Assay is used to measure the concentration of total GSH (reduced and oxidized) in a variety of animal and plant samples.

**Test Principle:**

The NWLSS™ Glutathione Assay is a modification of the method first described by Tietze.<sup>2</sup> The general thiol reagent, 5-5'-dithiobis[2-nitrobenzoic acid] (DTNB, Ellman's Reagent) reacts with GSH to form the 412 nm chromophore, 5-thionitrobenzoic acid (TNB) and GS-TNB. The GS-TNB is subsequently reduced by glutathione reductase and  $\beta$ -nicotinamide adenine dinucleotide phosphate (NADPH), releasing a second TNB molecule and recycling GSH; thus amplifying the response. Any oxidized GSH (GSSG) initially present in the reaction mixture or formed from the mixed disulfide reaction of GSH with GS-TNB is rapidly reduced to GSH.

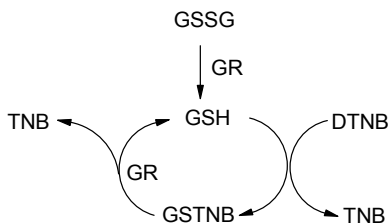

**Figure 2.**  
**GSH Recycling Pathway**

**General Specifications:**

Format:	96 wells	
Number of tests:	Triplicate =	29
	Duplicate =	45
Specificity:	Total Glutathione (GSH plus GSSG)	
Sensitivity:	LLD = 0.1 $\mu$ M in the reaction (Calculated as 3.29 standard deviations from zero).	

**Kit Contents:**

NADPH:	$\beta$ -Nicotinamide adenine dinucleotide phosphate, reduced; 1 bottle powder.	
NADPH Diluent:	Buffer with stabilizer; 1 x 6 mL vial.	
DTNB:	5-5'-Dithiobis(2-nitrobenzoic acid) in phosphate buffer with EDTA; 1 x 5.6 mL vial.	
GR Enzyme:	Glutathione reductase in Assay Buffer with protein stabilizer; 1 x 5.6 mL vial.	
Calibrators:	Glutathione disulfide in 10 mM HCl, 400 $\mu$ M GSH equivalents; 1 x 0.5 mL vial.	
Assay Buffer:	Phosphate buffer with EDTA, pH 7.6; 1 x 60 mL bottle.	

**Required Materials Not Provided:**

Pipettes capable of transferring 25, 50, 100 and 500  $\mu$ L volumes.  
A multi-channel or repeater pipette (recommended)  
Metaphosphoric acid, (Sigma 239275 or equivalent).  
Deionized water  
Polypropylene microcentrifuge tubes or equivalent  
Stir bar  
Beaker or flask 25-50 mL

**Required Instrumentation:**

Microplate reader with kinetics capability at 405 nm.

Note: The  $\lambda_{\text{max}}$  for TNB is 412 nm, however most plate readers are equipped with a 405 nm filter which will be fine for use with this assay

**Warnings, Limitations, Precautions:***NADPH*

$\beta$ -Nicotinamide adenine dinucleotide phosphate, reduced form (CAS 2646-71-1) is irritating to the eyes, respiratory system and skin. Target organs: nerves, liver. In case of contact with eyes, rinse immediately with copious amounts of water and seek medical advice. Wear suitable protective clothing.

*MPA*

Metaphosphoric acid (CAS 37267-86-0) is corrosive and may cause burns. In case of contact with eyes, rinse immediately with copious amounts of water and seek medical advice. Wear protective clothing, gloves and eye protection.

*DTNB*

5-5'-Dithiobis(2-nitrobenzoic acid) (CAS 69-78-3) is irritating to the eyes, respiratory system and skin. In case of contact with eyes, rinse immediately with copious amounts of water and seek medical advice. Wear suitable protective clothing.

*Interference*

Thiol containing compounds, such as cysteine,  $\beta$ -mercaptoethanol or dithiothreitol may compete with GSH for DTNB. N-ethylmaleimide or other thiol alkylating agents are known to interfere with GR and also will react with GSH and should be avoided.

**Storage Instructions:**

Upon receipt, store the reagents at 2-8°C. Do not use components beyond the expiration date printed on the label.

All reagents should be brought to room temperature (18-25°C) prior to use.

**Assay Preparation***Plate Reader Setup Recommendations*

Wavelength: 405 nm

Duration: 3 minutes

Interval: 15-20 seconds

Data Reduction: Linear regression

**Reagent Preparation:***5% Metaphosphoric acid*

1. Weigh 1 gram of MPA into a beaker or flask
2. Add 20 mL deionized water to the beaker or flask
3. Stir until dissolved
4. Store at 2-8 °C until use

Note: 5% MPA solution should be prepared and used the same day.

*NADPH*

1. Add entire contents of NADPH Diluent bottle to NADPH bottle.
2. Secure cap on bottle containing the reconstituted NADPH and mix briefly by inversion.

Note: Once reconstituted, the NADPH is usable for up to 1 week if stored at 2-8°C in the original container. The NADPH reagent can also be stored at -20 °C for up to 6 weeks. However, it is recommended that the reagent be tested before committing samples.

*Calibrators*

Add 50 µL GSH Calibrator to 950 µL Assay Buffer = 20 µM GSH Equiv.

Add 25 µL GSH Calibrator to 975 µL Assay Buffer = 10 µM GSH Equiv.

Assay Buffer Only = 0 µM GSH equivalents.

*Other Reagents*

The DTNB, GR and Assay Buffer are supplied ready-to-use.

**Sample Handling/Preparation:**

Note 1: If GSH is to be normalized to hemoglobin or protein, these tests must be performed on non-deproteinated samples.

Note 2: If there is a possibility that GSSG will be tested, customers should keep in mind that GSH oxidizes rapidly when in a non-reducing environment. Therefore, samples must be treated with a suitable GSH sequestering reagent in a manner and timeframe that seeks to minimize levels of artifact GSSG produced.

The multi-disciplinary interest in measuring GSH has resulted in a myriad of sample types and experimental conditions and is beyond the scope of this product insert to describe sample processing in detail. However, general guidelines are provided for representative sample types.

*Deproteinaton:*

Most samples will require removal of interfering proteins and metabolizing enzymes prior to assay. The acidic environment of the deproteinated sample also improves the stability of GSH.

**Sample Handling/Preparation (continued):***Deproteination (continued):*

1. Add 50  $\mu\text{L}$  of sample to microcentrifuge tube.
2. Add 100  $\mu\text{L}$  cold 5% MPA.
3. Vortex for 5 seconds (5 count).
4. Centrifuge at 10,000  $xg$  for 5 minutes or equivalent.
5. Carefully remove supernatant.
6. Place supernatant on ice for same day assay else store at  $-70^{\circ}\text{C}$ .
7. Just prior to assay, dilute sample at least 1/20 in Assay Buffer

*Whole Blood:*

Collect samples using EDTA, heparin, citrate or ACD anticoagulant. Samples should be kept cold ( $2-8^{\circ}\text{C}$ ) during transport and processing and should be deproteinated as soon as possible. If not assayed within a short time, sample extracts may be frozen at  $-70^{\circ}\text{C}$  for long term storage.

*Plasma:*

In most model systems plasma levels of GSH are typically at or below the level of detection of this assay. However, it is possible to modify the standard protocol to yield greater sensitivity. See [www.nwllifescience.com/products/assaykit/gshhighsens.htm](http://www.nwllifescience.com/products/assaykit/gshhighsens.htm) for additional information regarding the required modifications for measurement of plasma GSH.

*Tissue*

The GSH concentration in most tissues ranges from 1-10 mM; therefore, a 10% w/v homogenate is recommended. Contaminating blood contains significant concentrations of GSH and, if practical, should be removed by perfusion with an appropriate buffer, such as phosphate buffered saline containing heparin. Clarify homogenates by centrifugation and store the supernatant on ice. The GSH concentration can be normalized to the wet weight of the tissue sample or to the protein concentration of the homogenate. Deproteiniate as soon as possible to reduce the loss of GSH to various metabolic enzymes and mixed disulfide reactions. If not assayed within a short time, freeze both the homogenates and extracts at  $-70^{\circ}\text{C}$ .

*Cultured Cells*

Prepare a cell suspension of approximately  $10^6$  cells per mL in an appropriate buffer and homogenize or disrupt the cells by sonication. Clarify the homogenate by centrifugation and store the supernatant on ice. The GSH concentration can be normalized to the cell number or to the protein concentration of the homogenate. Deproteiniate as soon as possible to reduce the loss of GSH to various metabolic enzymes and mixed disulfide reactions. It is recommended that a trial assay with a representative sample be tested to determine if the samples are within the dynamic range of the assay. If the test sample is below the sensitivity of the standard procedure, please contact NWLSS for assistance. If not assayed within a short time, freeze both the homogenates and extracts at  $-70^{\circ}\text{C}$ .

**Sample Handling/Preparation (continued):*****Oxidized Glutathione***

Oxidized GSH (GSSG) can be estimated using the NWLSS GSH Assay by first treating the sample with a reagent suitable for sequestering reduced GSH. One such reagent is 4-vinylpyridine (4-VP, Sigma V-3877 or equivalent). The GSSG in most samples is at or below the detection limit of the standard assay method. However, a modified protocol for testing GSSG in tissue samples using 4-VP is available on our website at <http://www.nwlifescience.com/products/PDF/gssgapp07102013.htm>. N-ethylmaleimide (NEM) is another such sequestering reagent that can yield better results when used properly. For additional information on this technique, please contact our technical service department at [nwtech@nwlifescience.com](mailto:nwtech@nwlifescience.com)

**Assay Protocol:*****Standard Procedure***

1. Bring all reagents to room temperature.
2. Remove microplate from plastic bag.
3. Add 50  $\mu$ L of calibrator, diluted sample, and diluted control to a designated well.
4. Add 50  $\mu$ L of DTNB to each well.
5. Add 50  $\mu$ L of GR to each well
6. Incubate microplate for 3-5 minutes at room temperature.
7. Add 50  $\mu$ L of reconstituted NADPH to each well.
8. Begin recording the absorbance at 405 nm at 15-20 second intervals for 3 minutes.
9. Determine the concentration of the controls and samples. If using data reduction on the plate reader, skip steps a-d.
  - a. Calculate the rate for each calibrator, control and sample from the slope of the linear regression of  $A_{405}$  as a function of time.
  - b. Calculate the linear regression parameters to obtain the equation of the line.
  - c. Calculate the concentrations of the controls and samples.
  - d. Correct the control and sample for dilution and report results.

**Data Analysis:**


The following example shows the expected behavior and results of GSH determination using the standard method with the NWLSS™ Glutathione assay.

A 50  $\mu$ L sample of whole blood sample was added to a microcentrifuge tube containing 100  $\mu$ L cold 5% MPA. The microcentrifuge tube was then vortexed, centrifuged and the supernatant collected and placed on ice. Following a 25:500 dilution in Assay Buffer, the diluted sample was assayed using the standard method.



**Data Analysis (continued):**

The plot of the absorbance as a function of time for each calibrator and the unknown sample are shown in Figure 3. On inspection the curves are linear as expected.

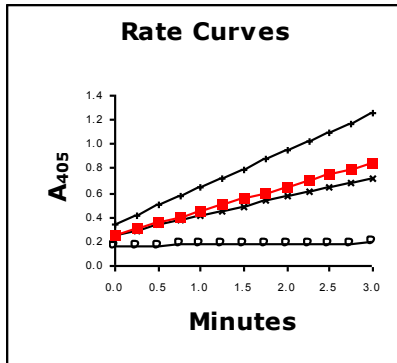

**Figure 3.**

Plot of the  $A_{405}$  as a function of time for calibrators 0  $\mu\text{M}$  ( $\circ$ ), 10  $\mu\text{M}$  ( $\times$ ), 20  $\mu\text{M}$  (+) GSH and unknown ( $\blacksquare$ ).

Alternatively, the  $r^2$  parameter can be used as a measure of linearity (Table 1). If curves appear to be non-linear or have an  $r^2$  value less than 0.990, then that sample should be repeated following further dilution, see Dynamic Range, page 11.

For each curve in Figure 3, the reaction rate (rate = slope) for each curve was determined using linear regression analysis. The rate for each calibrator and the unknown sample are shown in Table 1.

Sample		Rate $\Delta A_{405}/\text{min}$	$r^2$
Calibrator	0 $\mu\text{M}$	0.0088	0.996
Calibrator	10 $\mu\text{M}$	0.1560	1.000
Calibrator	20 $\mu\text{M}$	0.3016	1.000
Unknown		0.1972	1.000

**Table 1: Reaction Rates**

**Data Analysis: (continued)**

The calibration curve, shown in Figure 4, is constructed by plotting the rate for each calibrator as a function of the GSH concentration.


The general equation for the GSH calibration curve is:

$$\text{Rate} = a[\text{GSH}] + b \quad \text{Equation 1}$$

where  $a$  and  $b$  are the slope and intercept of the linear regression equation, respectively.

Linear regression analysis of the rates as a function of concentration from Table 1 yielded the following equation:

$$\text{Rate} = 0.01464[\text{GSH}] + 0.0091, r^2 = 1.000.$$

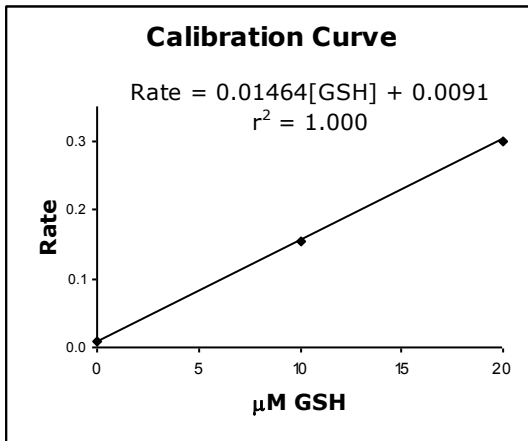

**Figure 4.**

Plot of the increase in the rate ( $\Delta A_{405}/\text{min}$ ) as a function of the concentration of GSH.

Rearranging Equation 1, the concentration of GSH in the unknown is determined given the rate...

$$[\text{GSH}] = \frac{\text{Rate} - b}{a} * df \quad \text{Equation 2}$$

where  $df$  is the dilution factor for the unknown.

**Data Analysis: (continued)**

Continuing with the example above...

- The unknown rate = 0.1972 A<sub>405</sub>/min
- The unknown *df* (1/3 deproteination and 1/21 buffer dilution) = 63x

Substituting the values into Equation 2, the GSH concentration for the unknown was found to be...

$$[\text{GSH}] = \frac{0.1972 - 0.0091}{0.0146} * 63 = 811 \mu\text{M}$$

**Performance Details:***Precision*

The precision of the assay was estimated by measuring a set of controls having low and high GSH concentrations, in duplicate, two times each day for five consecutive days using the standard procedure.

Control	μM GSH	Intra-assay		Inter-assay		Total	
	Mean	s	%CV	s	%CV	s	%CV
Low	219	6	2.7	9	4.1	7	3.2
High	827	11	1.3	14	1.7	12	1.5

*Sensitivity*


The method detection limit or MLD is estimated using twice the standard error of the rate of the calibration curve or 0.5 μM GSH. Alternatively, the lower limit of detection (LLD) is defined as 3.29 standard deviations from zero or 0.1 μM.

*Accuracy*

Recovery: A 10 mM solution of GSH was diluted in Assay Buffer to 20 μM and measured against the GSSG calibrators using the standard method. The GSH recovery was 102% showing the equivalency of the GSSG calibrators in the assay.

**Performance Details (continued):*****Dynamic Range:***

The linearity and sensitivity of the calibration curve defines the useful range of the assay. As shown in Figure 5, an unknown sample with a GSH concentration, following dilution in Assay Buffer, greater than the 20  $\mu\text{M}$  will not be accurate and should be further diluted and re-assayed.

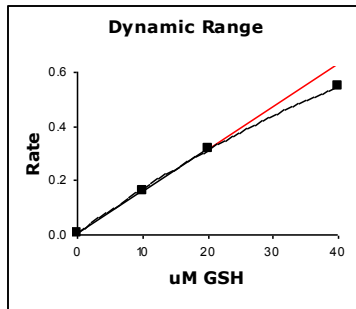

**Figure 5.**

The calibration curve (—) is compared to the plot of a series of GSH concentrations that exceed the range of the assay (■).

***Stability***

All unopened reagents are stable until the expiration date stated on the package label when stored at 2-8 °C.

Reconstituted NADPH is usable for up to 1 week if stored at 2-8 °C and up to 1 month at -20 °C if stored sealed and uncontaminated in the original bottle.

Deproteinized samples in 5% MPA are stable up to 6 months at -20 °C or indefinitely at -80 °C.

**Things to Note:**

1. Endpoint assay is not recommended; kinetic is more precise and faster.
2. The reagents are provided with a reasonable overfill but use caution if using multi-channel pipetting troughs to minimize reagent volume loss.
3. Do not pre-combine the NADPH, DTNB and GR. GR will catalyze the reduction of DTNB by NADPH resulting in higher background and/or exhaustion of reagents.

**Things to Note: (continued):**

4. NADPH that has been stored for a prolonged time can be tested by demonstrating that the rate curve for the high calibrator is linear over a 5 minute interval.
5. Unused wells can be protected from contamination by used wells by sealing the used wells with 1 inch cellophane tape.
6. GSH is rapidly metabolized and will undergo mixed disulfide reactions; therefore, it is recommended that samples be processed as soon as possible to avoid the loss of GSH.
7. Please contact NWLSS for advice if the expected GSH concentration in the sample is near or below the sensitivity of the standard method.

**References:**

1. Halliwell, B.; Gutteridge, J.M.C. Free Radicals in Biology and Medicine. Oxford University Press, New York; 1999.
2. Teitze, F. Enzymatic method for the quantitative determination of nanogram amounts of total and oxidized glutathione: Applications to mammalian blood and other tissues, *Anal. Biochem.***27**:502-522; 1969.


**Statement of Limited Warranty:**

Northwest Life Science Specialties, LLC (NWLSS) makes no guarantee of any kind, expressed or implied, that extends beyond the description of the material in this kit, except that they will meet our specifications at the time of delivery. Customer's remedy and NWLSS' sole liability is limited to, at NWLSS' option, refund of the purchase price, or the replacement of material not meeting our specification. By acceptance of our product, customer assumes all liability and will indemnify and hold NWLSS harmless for the consequence of this product's use or misuse by the customer, its employees, or others. Refund or replacement is conditioned of customer notifying NWLSS within twenty-one (21) days of the receipt of product. Failure to give notice within 21 days shall constitute a waiver by the customer of all claims hereunder with respect to said product.

**Procedure Checklist**

- \_\_\_ Process the sample and place on ice
- \_\_\_ Prepare 5% MPA and place on ice
- \_\_\_ Prepare Calibrators and place on ice
- \_\_\_ Bring GSH reagents to ambient temperature
- \_\_\_ Setup microplate reader
  - Wavelength = 405 nm
  - Mode = Kinetic
  - Reaction Time = 3 Minutes
  - Read Interval = 15-20 seconds (9-12 data points)
- \_\_\_ Deproteinase samples 50:100; sample to MPA
- \_\_\_ Vortex, 5 count
- \_\_\_ Centrifuge sample for 1 minute at 10,000 x g
- \_\_\_ Dilute samples 25:500; supernatant to Assay Buffer and place on ice
- \_\_\_ Reconstitute NADPH Reagent with NADPH Reagent Diluent
- \_\_\_ Add 50  $\mu$ L calibrators, controls and diluted samples to appropriate wells of microplate
- \_\_\_ Add 50  $\mu$ L DTNB Reagent to all wells
- \_\_\_ Add 50  $\mu$ L GR Reagent to all wells
- \_\_\_ Incubate 2-3 minutes at ambient temperature
- \_\_\_ Add 50  $\mu$ L NADPH Reagent
- \_\_\_ Place microplate in plate reader and begin 405 nm Measurements for 3 minutes
- \_\_\_ Calculate results
- \_\_\_ Return reagents to 2-8°C.

End-User Notes:

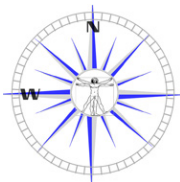

# **Northwest**

**Life Science Specialties, LLC**

**5131 NE 94th Avenue, Suite 201**

**Vancouver, WA 98662**

**Phone 360-449-3091 or Toll Free: 888-449-3091**

**Fax 360-449-3092**

**E-mail: [nwsales@nwlifescience.com](mailto:nwsales@nwlifescience.com)**